

## Product Datasheet

### Blood Group Lewis a Antibody [7LE] (orb1252414)

## Description

Blood Group Lewis a Antibody [7LE]

### Species/Host

Mouse

### Reactivity

Human, Mouse

### Conjugation

Unconjugated

### Tested

FC, IF, IHC-P, WB

### Applications

### Immunogen

Mucins isolated from ovarian cyst fluid were used as the immunogen for the Blood Group Lewis a antibody.

### Preservatives

PBS with 0.1 mg/ml rAlbumin and 0.05% sodium azide

### Form/Appearance

Liquid

### Concentration

0.2 mg/mL

### Storage

Aliquot and Store at 2-8°C. Avoid freeze-thaw cycles.

### Note

For research use only

### Application notes

Flow Cytometry: 0.5-1 ug/million cells in 0.1ml  
Immunofluorescence: 0.5-1 ug/ml  
Immunohistochemistry (FFPE): 0.5-1 ug/ml for 30 min at RT (1)Prediluted format: incubate for 30 min at RT (2)Optimal dilution of the Blood Group Lewis a antibody should be determined by the researcher.1. Staining of formalin-fixed tissues requires boiling tissue sections in 10mM Citrate buffer, pH 6.0, for 10-20 min followed by cooling at RT for 20 min.2. The prediluted format is supplied in a dropper bottle and is optimized for use in IHC. After epitope retrieval step (if required), drip mAb solution onto the tissue section and incubate at RT for 30 min.

### Isotype

IgG1, kappa

### Clonality


Monoclonal

### Uniprot ID


NotApplicable

### Dilution Range

Flow Cytometry: 0.5-1 ug/million cells in 0.1ml  
Immunofluorescence: 0.5-1 ug/ml  
Immunohistochemistry (FFPE): 0.5-1 ug/ml for 30 min at RT (1)Prediluted format: incubate for 30 min at RT (2)Optimal dilution of the Blood Group Lewis a antibody should be determined by the researcher.1. Staining of formalin-fixed tissues requires boiling tissue sections in 10mM Citrate buffer, pH 6.0, for 10-20 min followed by cooling at RT for 20 min.2. The prediluted format is supplied in a dropper bottle and is optimized for use in IHC. After epitope retrieval step (if required), drip mAb solution onto the tissue section and incubate

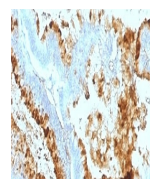

Formalin-fixed, paraffin-embedded human ...

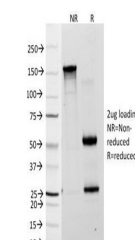

SDS-PAGE Analysis of Purified, BSA-Free ...