

Product Datasheet

TOX3 Antibody (orb1252546)

Description

TOX3 Antibody

Species/Host

Mouse

Reactivity

Human

Conjugation

Unconjugated

Tested

FC, IF, IHC-P, WB

Applications
Immunogen

A recombinant fragment (139 amino acid residues around aa 200-400) from the human protein was used as the immunogen for the TOX3 antibody.

Target

TOX3

Preservatives

PBS with 0.1 mg/ml rAlbumin and 0.05% sodium azide

Form/Appearance

Liquid

Concentration

0.2 mg/mL

Storage

Aliquot and Store at 2-8°C. Avoid freeze-thaw cycles.

Note

For research use only

Application notes

Flow Cytometry: 0.5-1 ug/million cells in 0.1ml
Immunofluorescence: 0.5-1 ug/ml
Immunohistochemistry (FFPE): 0.5-1 ug/ml for 30 min at RT (1)Prediluted format : incubate for 30 min at RT (2)Optimal dilution of the TOX3 antibody should be determined by the researcher.1. Staining of formalin-fixed tissues requires boiling tissue sections in 10mM Tris with 1mM EDTA, pH 9.0 or 10mM Citrate buffer, pH 6.0, for 10-20 min followed by cooling at RT for 20 min2. The prediluted format is supplied in a dropper bottle and is optimized for use in IHC. After epitope retrieval step (if required), drip mAb solution onto the tissue section and incubate at RT for 30 min.

Isotype

IgG2b, kappa

Clonality


Monoclonal

Uniprot ID


O15405

Dilution Range


Flow Cytometry: 0.5-1 ug/million cells in 0.1ml
Immunofluorescence: 0.5-1 ug/ml
Immunohistochemistry (FFPE): 0.5-1 ug/ml for 30 min at RT (1)Prediluted format : incubate for 30 min at RT (2)Optimal dilution of the TOX3 antibody should be determined by the researcher.1. Staining of formalin-fixed tissues requires boiling tissue sections in 10mM Tris with 1mM EDTA, pH 9.0 or 10mM Citrate buffer, pH 6.0, for 10-20 min followed by cooling at RT for 20 min2. The prediluted

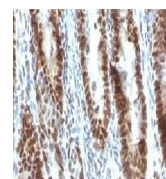

Formalin-fixed, paraffin-embedded human ...

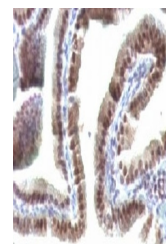

Formalin-fixed, paraffin-embedded human ...

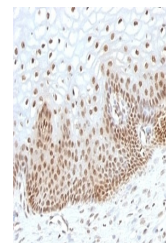

Formalin-fixed, paraffin-embedded human ...