

---

## Product Datasheet

### Nucleoli Marker Antibody [NM95] (orb1252651)

**Description**

Nucleoli Marker Antibody [NM95]

**Species/Host**

Mouse

**Reactivity**

Human

**Conjugation**

Unconjugated

**Tested Applications**

FC, ICC, IF, IHC

**Immunogen**

Nuclei of myeloid leukemia biopsy cells were used as the immunogen.

**Preservatives**

PBS with 0.1 mg/ml rAlbumin and 0.05% sodium azide

**Form/Appearance**

Liquid

**Concentration**

0.2 mg/mL

**Storage**

Aliquot and Store at 2-8°C. Avoid freeze-thaw cycles.

**Note**

For research use only

**Application notes**

Flow Cytometry: 0.5-1ug/10<sup>6</sup> cells  
 IF: 0.5-1 ug/ml  
 Immunocytochemistry (Acetone-fixed): 0.5-1 ug/ml for 30 min at RT  
 IHC (Formalin-paraffin): 0.5-1 ug/ml for 30 min at RT (1)  
 Prediluted format : incubate for 30 min at RT (2)  
 The concentration stated for each application is a general starting point. Variations in protocols, secondaries and substrates may require the Nucleoli marker antibody to be titrated up or down for optimal performance.  
 1. Staining of formalin-fixed tissues is enhanced by boiling tissue sections in 10mM Citrate Buffer, pH 6.0, for 10-20 min followed by cooling at RT for 20 minutes.  
 2. The prediluted format is supplied in a dropper bottle and is optimized for use in IHC. After epitope retrieval step (if required), drip mAb solution onto the tissue section and incubate at RT for 30 min.

**Isotype**


IgG1, kappa

**Clonality**


Monoclonal

**Dilution Range**


Flow Cytometry: 0.5-1ug/10<sup>6</sup> cells  
 IF: 0.5-1 ug/ml  
 Immunocytochemistry (Acetone-fixed): 0.5-1 ug/ml for 30 min at RT  
 IHC (Formalin-paraffin): 0.5-1 ug/ml for 30 min at RT (1)  
 Prediluted format : incubate for 30 min at RT (2)  
 The concentration stated for each application is a general starting point. Variations in protocols, secondaries and substrates may require the Nucleoli marker antibody to be titrated up or down for optimal

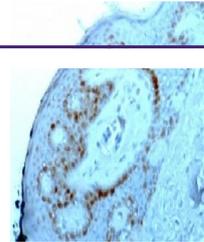

IHC staining of FFPE human skin with Nuc...

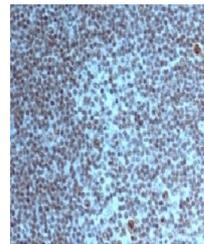

IHC staining of FFPE human tonsil tissue...

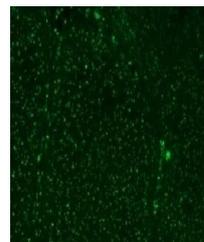

Immunofluorescent testing of colon carci...