

Product Datasheet

Golgi Complex Antibody [371-4] (orb1252662)

Description

Golgi Complex Antibody [371-4]

Species/Host

Mouse

Reactivity

Human

Conjugation

Unconjugated

Tested

FC, ICC, IF, IHC, WB

Applications

Immunogen

SU-DHL-1 large cell lymphoma cells were used as the immunogen.

Preservatives

PBS with 0.1 mg/ml rAlbumin and 0.05% sodium azide

Form/Appearance

Liquid

Concentration

0.2 mg/mL

Storage

Aliquot and Store at 2-8°C. Avoid freeze-thaw cycles.

Note

For research use only

Application notes

Western Blot: 0.5-1 ug/ml Flow Cytometry: 0.5-1 ug/million cells in 0.1ml Immunofluorescence: 0.5-1 ug/ml Immunocytochemistry: (Acetone/paraformaldehyde fixed) 0.5-1 ug/ml for 30 min IHC (FFPE): 0.5-1 ug/ml for 30 minutes at RT (1) Prediluted format : incubate for 30 min at RT (2) The concentration stated for each application is a general starting point. Variations in protocols, secondaries and substrates may require the Golgi complex antibody to be titrated up or down for optimal performance. 1. Staining of formalin-fixed tissues requires boiling tissue sections in 10mM citrate buffer, pH 6.0, for 10-20 min followed by cooling at RT for 20 minutes. 2. The prediluted format is supplied in a dropper bottle and is optimized for use in IHC. After epitope retrieval step (if required), drip mAb solution onto the tissue section and incubate at RT for 30 min.

Isotype


IgG1, kappa

Clonality

Monoclonal

Dilution Range

Western Blot: 0.5-1 ug/ml Flow Cytometry: 0.5-1 ug/million cells in 0.1ml Immunofluorescence: 0.5-1 ug/ml Immunocytochemistry: (Acetone/paraformaldehyde fixed) 0.5-1 ug/ml for 30 min IHC (FFPE): 0.5-1 ug/ml for 30 minutes at RT (1) Prediluted format : incubate for 30 min at RT (2) The concentration stated for each application is a general starting point. Variations in protocols, secondaries and substrates may require the Golgi complex antibody to be titrated up or down for optimal performance. 1

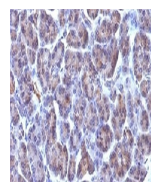

Formalin-fixed, paraffin-embedded human ...



SDS-PAGE Analysis of Purified, BSA-Free ...