

Product Datasheet

General Glucagon Like Peptide 1 (GLP1) ELISA Kit (orb1088181)

Catalog Number orb1088181

Description This assay employs the competitive inhibition enzyme immunoassay technique.

The microtiter plate provided in this kit has been pre-coated with Glucagon Like Peptide 1(GLP1) protein. Standards or samples are added to the appropriate microtiter plate wells then with a biotin-conjugated antibody specific to Glucagon Like Peptide 1(GLP1). Next, Avidin conjugated to Horseradish Peroxidase (HRP) is added to each microplate well and incubated. After TMB substrate solution is added. The enzyme-substrate reaction is terminated by the addition of sulphuric acid solution and the color change is measured spectrophotometrically at a wavelength of $450 \, \text{nm} \pm 10 \, \text{nm}$. The concentration of Glucagon Like Peptide 1(GLP1) in the samples is then determined by comparing the OD of the samples

to the standard curve.

Reactivity All

Range 15.63-1000 pg/mL

Concentration 1000 pg/mL

Note For research use only

Application notes standard: 1000 pg/mL. Test principle: This assay employs the competitive

kit has been pre-coated with GLP1 protein. Standards or samples are added to the appropriate microtiter plate wells then with a biotin-conjugated antibody specific to GLP1. Next, Avidin conjugated to Horseradish Peroxidase (HRP) is added to each microplate well and incubated. After TMB substrate solution is added. The enzyme-substrate reaction is terminated by the addition of sulphuric acid solution and the color change is measured spectrophotometrically at a

wavelength of 450nm \pm 10nm. The concentration of GLP1 in the samples is then

inhibition enzyme immunoassay technique. The microtiter plate provided in this

determined by comparing the OD of the samples to the standard curve

Sample Types serum, plasma and other biological fluids






Assay Time 2.5h

Sensitivity 4.18 pg/mL

Expiration Date Please enquire.

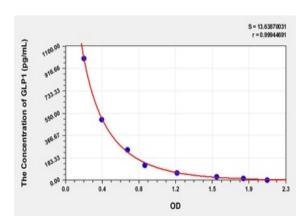

Email: <u>info@biorbyt.com</u>, <u>support@biorbyt.com</u>
Phone: <u>+1 (415) 906-5211</u> | Fax: <u>+1 (415) 651-8558</u>