

Product Datasheet

Anti-Cyclin B1 Antibody (orb1423587)

Description Rabbit polyclonal antibody to Cyclin B1.

Species/Host Rabbit

Reactivity Human

Conjugation Unconjugated

Tested Applications IF, IHC, WB

Storage Maintain refrigerated at 2-8°C for up to 2 weeks. For long term storage store at -

20°C in small aliquots to prevent freeze-thaw cycles

Note For research use only

Clonality Polyclonal


Clone Number CCNB1

Uniprot ID P14635

Dilution Range WB: WB (1/500 - 1/1000), IH (1/100 - 1/200), IF/IC (1/100 - 1/500), IF: WB (1/500 -

1/1000), IH (1/100 - 1/200), IF/IC (1/100 - 1/500)

Expiration Date 12 months from date of receipt.

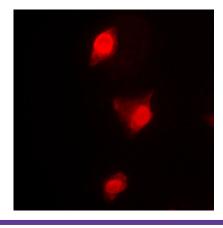

Immunofluorescent analysis of Cyclin B1 staining in K562 Cells. Formalin-fixed cells were permeabilized with 0.1% Triton X-100 in TBS for 5-10 minutes and blocked with 3% BSA-PBS for 30 minutes at room temperature. Cells were probed with the primary antibody in 3% BSA-PBS and incubated overnight at 4°C in a humidified chamber. Cells were washed with PBST and incubated with a DyLight 594-conjugated secondary antibody (red) in PBS at room temperature in the dark.

Biorbyt Ltd.


7 Signet Court, Swann's Road, Cambridge, CB5 8LA, United Kingdom Email: info@biorbyt.com Phone: +44 (0) 1223 859-353 | Fax: +1 (415) 651-8558

Biorbyt LLC.


68 TW Alexander Drive,
Durham, NC, 27713, United States
Email: info@biorbyt.com
Phone: +1 (415) 906-5211 | Fax: +1 (415) 651-8558





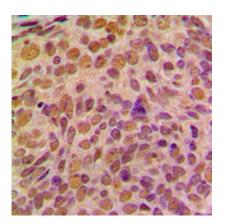

Immunohistochemical analysis of Cyclin B1 staining in human breast cancer formalin fixed paraffin embedded tissue section. The section was pre-treated using heat mediated antigen retrieval with sodium citrate buffer (pH 6.0). The section was then incubated with the antibody at room temperature and detected using an HRP conjugated compact polymer system. DAB was used as the chromogen. The section was then counterstained with haematoxylin and mounted with DPX.

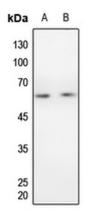

Western blot analysis of Cyclin B1 expression in HEK293T (A), Hela (B) whole cell lysates.

Phone: <u>+1 (415) 906-5211</u> | Fax: <u>+1 (415) 651-8558</u>