

Product Datasheet

TUBB3 (Neuronal Marker) Rabbit Polyclonal Antibody (BF594) (orb1587539)


Description	TUBB3 (Neuronal Marker) Rabbit Polyclonal Antibody (BF594)
Species/Host	Rabbit
Reactivity	Human, Mouse, Rat
Conjugation	BF594
Tested Applications	FC, ICC, IF
Immunogen	KLH conjugated synthetic peptide derived from human beta III Tubulin (401-450/450aa)
Target	TUBB3
Preservatives	0.01M TBS (pH7.4) with 1% rAlbumin, 0.02% Proclin300 and 50% Glycerol.
Concentration	1mg/ml
Storage	Shipped at 4°C. Store at -20°C for one year. Avoid repeated freeze/thaw cycles.
Note	For research use only
Isotype	IgG
Clonality	Polyclonal
Antibody Type	Primary Antibody
Dilution Range	Flow-Cyt=1µg/Test, IF=1:200-800, ICC/IF=1:100
Expiration Date	12 months from date of receipt.

Biorbyt Ltd.


7 Signet Court, Swann's Road,
Cambridge, CB5 8LA, United Kingdom
Email: info@biorbyt.com, support@biorbyt.com
Phone: [+44 \(0\) 1223 859-353](tel:+441223859353) | Fax: [+1 \(415\) 651-8558](tel:+14156518558)

Biorbyt LLC.

68 TW Alexander Drive,
Durham, NC, 27713, United States
Email: info@biorbyt.com, support@biorbyt.com
Phone: [+1 \(415\) 906-5211](tel:+14159065211) | Fax: [+1 \(415\) 651-8558](tel:+14156518558)

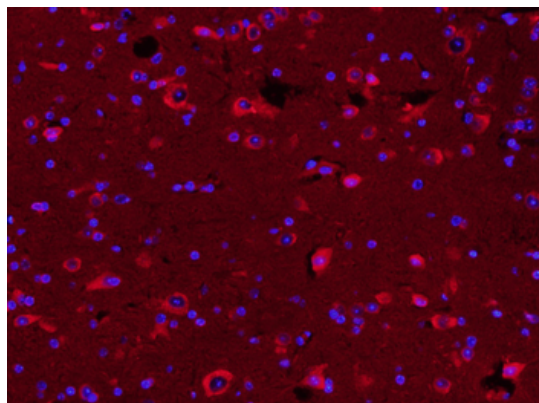

Paraformaldehyde-fixed, paraffin embedded (Mouse brain), Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15 min, Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes, Blocking buffer (normal goat serum) at 37°C for 30 min, Incubation: Anti-TUBB3 Antibody, conjugated 1:400, 1.5 hours at 37°C, DAPI (5 ug/ml, blue) was used to stain the cell nuclei.

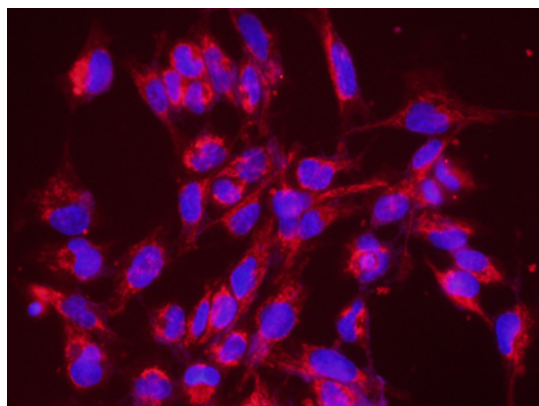

Tissue/Cell: U251 cell, 4% Paraformaldehyde-fixed, Triton X-100 at room temperature for 20 min, Blocking buffer (normal goat serum) at 37°C for 20 min, Incubation: Anti-TUBB3 (Neuronal Marker) Polyclonal Antibody, conjugated (orb1587539) 1:200, 1.5 hours at 37°C, DAPI (5 ug/ml, blue) was used to stain the cell nuclei.

Biorbyt Ltd.

7 Signet Court, Swann's Road,
Cambridge, CB5 8LA, United Kingdom
Email: info@biorbyt.com, support@biorbyt.com
Phone: [+44 \(0\) 1223 859-353](tel:+44(0)1223859353) | Fax: [+1 \(415\) 651-8558](tel:+1(415)6518558)

Biorbyt LLC.

68 TW Alexander Drive,
Durham, NC, 27713, United States
Email: info@biorbyt.com, support@biorbyt.com
Phone: [+1 \(415\) 906-5211](tel:+1(415)9065211) | Fax: [+1 \(415\) 651-8558](tel:+1(415)6518558)