

## Product Datasheet

### Perforin Antibody (APC) (orb1742786)

<b>Catalog Number</b>	orb1742786
<b>Category</b>	Antibodies
<b>Description</b>	Perforin Antibody (APC)
<b>Target</b>	Perforin
<b>Clonality</b>	Monoclonal
<b>Isotype</b>	Mouse IgG2b kappa
<b>Conjugation</b>	APC
<b>Reactivity</b>	Bovine, Human
<b>Buffer/Preservatives</b>	Stabilizing phosphate buffered saline (PBS), pH 7.4, 15 mM sodium azide
<b>Purification</b>	Purified antibody is conjugated with activated allophycocyanin (APC) under optimum conditions and unconjugated antibody and free fluorochrome are removed by size-exclusion chromatography.
<b>Immunogen</b>	purified granules from human YT lymphoma cell line
<b>UniProt ID</b>	<b>P14222</b>
<b>Tested applications</b>	FC
<b>Application notes</b>	Flow cytometry: The reagent is designed for analysis of human blood cells using 10 µl reagent / 100 µl of whole blood or 10 <sup>6</sup> cells in a suspension. The content of a vial (1 ml) is sufficient for 100 tests. Intracellular staining.
<b>Specificity</b>	The mouse monoclonal antibody dG9 (also known as deltaG9) recognizes perforin, a 70 kDa protein expressed in cytoplasmic granules of cytotoxic T cells and NK cells.

**Biorbyt Ltd.**

7 Signet Court, Swann Road  
Cambridge  
CB5 8LA  
United Kingdom


Email: [info@biorbyt.com](mailto:info@biorbyt.com), [support@biorbyt.com](mailto:support@biorbyt.com)  
Phone: [+44 \(0\)1223 859353](tel:+44(0)1223859353) | Fax: [+1 \(415\) 651-8558](tel:+1(415)651-8558)

**Biorbyt LLC**


68 TW Alexander Drive  
Research Triangle Park  
Durham  
NC 27713-2847  
United States

Email: [info@biorbyt.com](mailto:info@biorbyt.com), [support@biorbyt.com](mailto:support@biorbyt.com)  
Phone: [+1 \(415\) 906-5211](tel:+1(415)906-5211) | Fax: [+1 \(415\) 651-8558](tel:+1(415)651-8558)

<b>Clone Number</b>	dG9
<b>Storage</b>	Store at 2-8°C. Protect from prolonged exposure to light. Do not freeze.
<b>Note</b>	For research use only
<b>Entrez</b>	<b>5551</b>
<b>Expiration Date</b>	12 months from date of receipt.

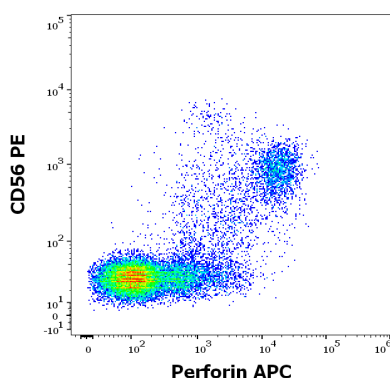

Flow cytometry multicolor surface staining pattern of human lymphocytes using anti-human CD56 (LT56) PE antibody (10 µl reagent / 100 µl of peripheral whole blood) and intracellular staining using anti-Perforin (dG9) APC antibody (10 µl reagent / 100 µl of peripheral whole blood).

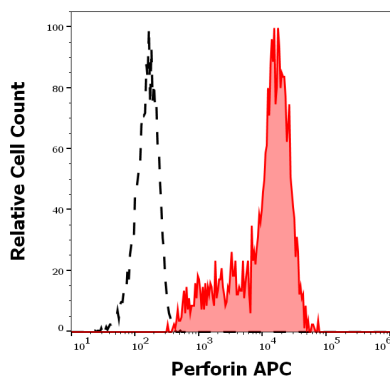

Separation of human Perforin positive CD56 positive lymphocytes (red-filled) from neutrophil granulocytes (black-dashed) in flow cytometry analysis (intracellular staining) of human peripheral whole blood stained using anti-Perforin (dG9) APC antibody (10 µl reagent / 100 µl of peripheral whole blood).

**Biorbyt Ltd.**

7 Signet Court, Swann Road  
Cambridge  
CB5 8LA  
United Kingdom

Email: [info@biorbyt.com](mailto:info@biorbyt.com), [support@biorbyt.com](mailto:support@biorbyt.com)


Phone: +44 (0)1223 859353 | Fax: +1 (415) 651-8558

**Biorbyt LLC**

68 TW Alexander Drive  
Research Triangle Park  
Durham  
NC 27713-2847  
United States

Email: [info@biorbyt.com](mailto:info@biorbyt.com), [support@biorbyt.com](mailto:support@biorbyt.com)

Phone: +1 (415) 906-5211 | Fax: +1 (415) 651-8558

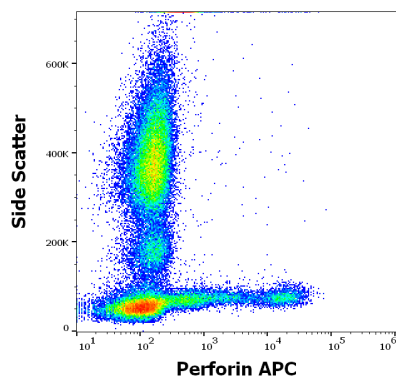

Flow cytometry intracellular staining pattern of human peripheral whole blood stained using anti-Perforin (dG9) APC antibody (10 µl reagent / 100 µl of peripheral whole blood).

**Biorbyt Ltd.**

7 Signet Court, Swann Road  
Cambridge  
CB5 8LA  
United Kingdom

Email: [info@biorbyt.com](mailto:info@biorbyt.com), [support@biorbyt.com](mailto:support@biorbyt.com)

Phone: [+44 \(0\)1223 859353](tel:+44(0)1223859353) | Fax: [+1 \(415\) 651-8558](tel:+1(415)651-8558)

**Biorbyt LLC**

68 TW Alexander Drive  
Research Triangle Park  
Durham  
NC 27713-2847  
United States

Email: [info@biorbyt.com](mailto:info@biorbyt.com), [support@biorbyt.com](mailto:support@biorbyt.com)

Phone: [+1 \(415\) 906-5211](tel:+1(415)906-5211) | Fax: [+1 \(415\) 651-8558](tel:+1(415)651-8558)