

Product Datasheet

VAMP4 Antibody (Center) (orb1930320)

Description	Affinity Purified Rabbit Polyclonal Antibody (Pab)
Species/Host	Rabbit
Reactivity	Human, Mouse
Conjugation	Unconjugated
Tested Applications	IHC-P, WB
Target	This VAMP4 antibody is generated from rabbits immunized with a KLH conjugated synthetic peptide between 16-45 amino acids from the Central region of human VAMP4.
Form/Appearance	Purified polyclonal antibody supplied in PBS with 0.09% (W/V) sodium azide. This antibody is purified through a protein A column, followed by peptide affinity purification.
Storage	Maintain refrigerated at 2-8°C for up to 2 weeks. For long term storage store at -20°C in small aliquots to prevent freeze-thaw cycles
Note	For research use only
Isotype	Rabbit IgG
Clonality	Polyclonal
Clone Number	RB22111
Antibody Type	Primary Antibody
MW	16397 Da
Uniprot ID	O75379
NCBI	NP_001172056.1, NP_003753.2

Biorbyt Ltd.

7 Signet Court, Swann's Road,
Cambridge, CB5 8LA, United Kingdom
Email: info@biorbyt.com, support@biorbyt.com
Phone: [+44 \(0\) 1223 859-353](tel:+44201223859353) | Fax: [+1 \(415\) 651-8558](tel:+14156518558)

Biorbyt LLC.


68 TW Alexander Drive,
Durham, NC, 27713, United States
Email: info@biorbyt.com, support@biorbyt.com
Phone: [+1 \(415\) 906-5211](tel:+14159065211) | Fax: [+1 \(415\) 651-8558](tel:+14156518558)

Dilution Range


WB: 1:1000, WB: 1:1000, IHC-P: 1:10~50

Expiration Date


12 months from date of receipt.

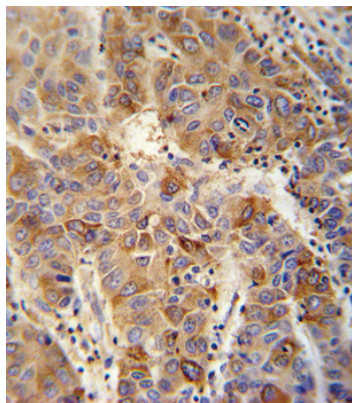

VAMP4 antibody (Center) immunohistochemistry analysis in formalin fixed and paraffin embedded human hepatocarcinoma followed by peroxidase conjugation of the secondary antibody and DAB staining. This data demonstrates the use of the VAMP4 antibody (Center) for immunohistochemistry. Clinical relevance has not been evaluated.

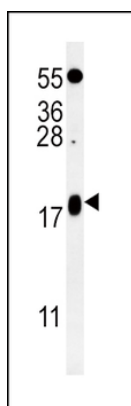

VAMP4 Antibody (Center) western blot analysis in K562 cell line lysates (35 ug/lane). This demonstrates the VAMP4 antibody detected the VAMP4 protein (arrow).

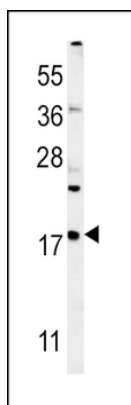

VAMP4 Antibody (Center) western blot analysis in mouse liver tissue lysates (35 ug/lane). This demonstrates the VAMP4 antibody detected the VAMP4 protein (arrow).

Biorbyt Ltd.

7 Signet Court, Swann's Road,
Cambridge, CB5 8LA, United Kingdom
Email: info@biorbyt.com, support@biorbyt.com
Phone: [+44 \(0\) 1223 859-353](tel:+44201223859353) | Fax: [+1 \(415\) 651-8558](tel:+14156518558)

Biorbyt LLC.

68 TW Alexander Drive,
Durham, NC, 27713, United States
Email: info@biorbyt.com, support@biorbyt.com
Phone: [+1 \(415\) 906-5211](tel:+14159065211) | Fax: [+1 \(415\) 651-8558](tel:+14156518558)