

## **Product Datasheet**

## Caldesmon / CAD Antibody HMW (orb322968)

Catalog Number orb322968

**Category** Antibodies

**Description** Caldesmon is an actin- and myosin-binding protein implicated in the regulation

of actomyosin interactions in smooth muscle and nonmuscle cells (could act as a

bridge between myosin and actin filaments). Stimulates actin binding of tropomyosin which increases the stabilization of actin filament structure. In muscle tissues, inhibits the actomyosin ATPase by binding to F-actin. This

inhibition is attenuated by calcium-calmodulin and is potentiated by

tropomyosin. Interacts with actin, myosin, two molecules of tropomyosin and with calmodulin. Also play an essential role during cellular mitosis and receptor

capping. Involved in Schwann cell migration during peripheral nerve

regeneration (By similarity). [UniProt]

**Clonality** Recombinant

**Species/Host** Rabbit

**Isotype** Rabbit IgG, kappa

**Conjugation** Unconjugated

**Reactivity** Human

**Immunogen** A full length human protein was used as the immunogen for the recombinant

CAD antibody.

UniProt ID Q05682

**Tested applications** IHC-P

**Dilution range** Immunohistochemistry (FFPE): 1-2ug/ml for 30 min at RT, Prediluted IHC only

format: incubate for 30 min at RT (1)





**Application notes**Titering of the recombinant CAD antibody may be required for optimal

performance.1. The prediluted format is supplied in a dropper bottle and is optimized for use in IHC. After epitope retrieval step (if required), drip mAb

solution onto the tissue section and incubate at RT for 30 min.

**Antibody Type** Primary Antibody


Clone Number RMCDN1-1

Formula 0.2 mg/ml in 1X PBS with 0.1 mg/ml BSA (US sourced) and 0.05% sodium azide

**Storage** Maintain refrigerated at 2-8°C for up to 2 weeks. For long term storage store at -

20°C in small aliquots to prevent freeze-thaw cycles.

**Note** For research use only

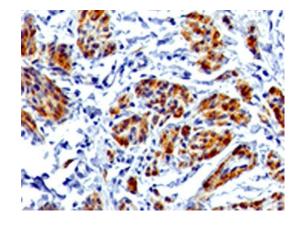

IHC testing of FFPE human uterus with recombinant CAD antibody (clone RMCDN1-1). HIER: steam sections in 10 mM Tris with 1 mM EDTA, pH9, for 10-20 min.

Phone: +1 (415) 906-5211 | Fax: +1 (415) 651-8558