



Product Datasheet Anti-TBX4 Antibody (orb745928)

Catalog Number	orb745928
Description	Anti-TBX4 Antibody. Tested in Flow Cytometry, IF, ICC, WB applications. This antibody reacts with Human.
Species/Host	Rabbit
Reactivity	Human
Conjugation	Unconjugated
Tested Applications	FC, ICC, IF, WB
Immunogen	A synthetic peptide corresponding to a sequence in the middle region of human TBX4, which shares 81.8% amino acid (aa) sequence identity with mouse TBX4.
Form/Appearance	Lyophilized
Concentration	Adding 0.2 ml of distilled water will yield a concentration of 500 μ g/ml.
Storage	Maintain refrigerated at 2-8°C for up to 2 weeks. For long term storage store at - 20°C in small aliquots to prevent freeze-thaw cycles.
Note	For research use only
Application notes	Western blot, 0.25-0.5µg/ml, Human Immunocytochemistry/Immunofluorescence, 5µg/ml, Human Flow Cytometry (Fixed), 1-3µg/1x106 cells, Human. Add 0.2ml of distilled water will yield a concentration of 500ug/ml
lsotype	Rabbit IgG
Clonality	Polyclonal
Antibody Type	Primary Antibody
MW	60 kDa

Biorbyt Ltd.

7 Signet Court, Swann's Road, Cambridge, CB5 8LA, United Kingdom Email: <u>info@biorbyt.com</u>, <u>support@biorbyt.com</u> Phone: <u>+44 (0) 1223 859-353</u> | Fax: <u>+1 (415) 651-8558</u>

Biorbyt LLC.

68 TW Alexander Drive, Durham, NC, 27713, United States Email: <u>info@biorbyt.com</u>, <u>support@biorbyt.com</u> Phone: <u>+1 (415) 906-5211</u> | Fax: <u>+1 (415) 651-8558</u>




Biorbyt.com

Uniprot ID


P57082

Expiration Date

12 months from date of receipt.

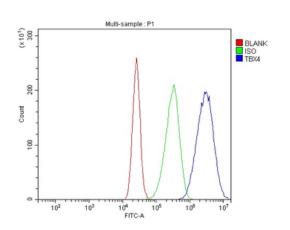

Flow Cytometry analysis of CACO-2 cells using anti-TBX4 antibody. Overlay histogram showing CACO-2 cells (Blue line). To facilitate intracellular staining, cells were fixed with 4% paraformaldehyde and permeabilized with permeabilization buffer. The cells were blocked with 10% normal goat serum. And then incubated with rabbit anti-TBX4 Antibody (1 μ g/1x10^6 cells) for 30 min at 20°C. DyLight®488 conjugated goat anti-rabbit IgG (5-10 μ g/1x10^6 cells) was used as secondary antibody for 30 minutes at 20°C. Isotype control antibody (Green line) was rabbit IgG (1 μ g/1x10^6) used under the same conditions. Unlabelled sample without incubation with primary antibody and secondary antibody (Red line) was used as a blank control.

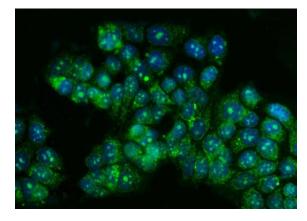

IF analysis of TBX4 using anti-TBX4 antibody. TBX4 was detected in immunocytochemical section of MCF-7 cells. Enzyme antigen retrieval was performed using IHC enzyme antigen retrieval reagent for 15 mins. The cells were blocked with 10% goat serum. And then incubated with 5 μg/mL rabbit anti-TBX4 Antibody overnight at 4°C. DyLight®488 Conjugated Goat Anti-Rabbit IgG was used as secondary antibody at 1:100 dilution and incubated for 30 minutes at 37°C. The section was counterstained with DAPI. Visualize using a fluorescence microscope and filter sets appropriate for the label used.

Biorbyt Ltd.


7 Signet Court, Swann's Road, Cambridge, CB5 8LA, United Kingdom Email: <u>info@biorbyt.com</u>, <u>support@biorbyt.com</u> Phone: +44 (0) 1223 859-353 | Fax: +1 (415) 651-8558

Biorbyt LLC.

68 TW Alexander Drive, Durham, NC, 27713, United States Email: <u>info@biorbyt.com</u>, <u>support@biorbyt.com</u> Phone: <u>+1 (415) 906-5211</u> | Fax: <u>+1 (415) 651-8558</u>



Biorbyt.com

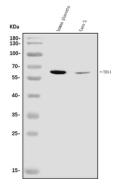

Western blot analysis of TBX4 using anti-TBX4 antibody. Electrophoresis was performed on a 5-20% SDS-PAGE gel at 70V (Stacking gel) / 90V (Resolving gel) for 2-3 hours. The sample well of each lane was loaded with 30 ug of sample under reducing conditions. Lane 1: human placenta tissue lysates, Lane 2: human Caco-2 whole cell lysates. After Electrophoresis, proteins were transferred to a Nitrocellulose membrane at 150 mA for 50-90 minutes. Blocked the membrane with 5% Non-fat Milk/ TBS for 1.5 hour at RT. The membrane was incubated with rabbit anti-TBX4 antigen affinity purified polyclonal antibody at 0.5 µg/mL overnight at 4°C, then washed with TBS-0.1% Tween 3 times with 5 minutes each and probed with a goat anti-rabbit IgG-HRP secondary antibody at a dilution of 1:5000 for 1.5 hour at RT. The signal is developed using an Enhanced Chemiluminescent detection (ECL) kit with Tanon 5200 system. A specific band was detected for TBX4 at approximately 60 KD. The expected band size for TBX4 is at 60 KD.

Biorbyt Ltd.

7 Signet Court, Swann's Road, Cambridge, CB5 8LA, United Kingdom Email: <u>info@biorbyt.com</u>, <u>support@biorbyt.com</u> Phone: +44 (0) 1223 859-353 | Fax: +1 (415) 651-8558

Biorbyt LLC.

68 TW Alexander Drive, Durham, NC, 27713, United States Email: <u>info@biorbyt.com</u>, <u>support@biorbyt.com</u> Phone: <u>+1 (415) 906-5211</u> | Fax: <u>+1 (415) 651-8558</u>