

Product Datasheet

Anti-Ki67 Antibody (monoclonal, 5E12) (orb763091)

Description Anti-Ki67 Antibody (monoclonal, 5E12). Tested in Flow Cytometry, IF, IHC, ICC

applications. This antibody reacts with Human.

Species/Host Mouse

Reactivity Human

Conjugation Unconjugated

Tested Applications FC, ICC, IF, IHC

Immunogen E. coli-derived human Ki67 recombinant protein (Position: K2860-I3256).

Form/Appearance Lyophilized

Concentration Adding 0.2 ml of distilled water will yield a concentration of 500 μg/ml.

Storage Maintain refrigerated at 2-8°C for up to 2 weeks. For long term storage store at -

20°C in small aliquots to prevent freeze-thaw cycles.

Note For research use only

Application notes Immunohistochemistry (Paraffin-embedded Section), 2-5µg/ml, Human

Immunocytochemistry/Immunofluorescence, 5µg/ml, Human

Immunofluorescence, 5µg/ml, Human Flow Cytometry (Fixed), 1-3µg/1x106 cells,

Human. Add 0.2ml of distilled water will yield a concentration of 500ug/ml

Isotype Mouse IgG2b

Clonality Monoclonal

Clone Number 5E12

Antibody Type Primary Antibody







MW 35-60 kDa

Uniprot ID P46013

Expiration Date 12 months from date of receipt.

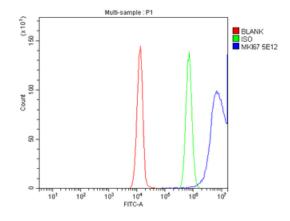

Flow Cytometry analysis of Jurkat cells using anti-Ki67 antibody. Overlay histogram showing Jurkat cells (Blue line). The cells were blocked with 10% normal goat serum. And then incubated with mouse anti-Ki67 Antibody (1 μ g/1x10^6 cells) for 30 min at 20°C. DyLight® 488 conjugated goat anti-mouse IgG (5-10 μ g/1x10^6 cells) was used as secondary antibody for 30 minutes at 20°C. Isotype control antibody (Green line) was mouse IgG (1 μ g/1x10^6) used under the same conditions. Unlabelled sample (Red line) was also used as a control.

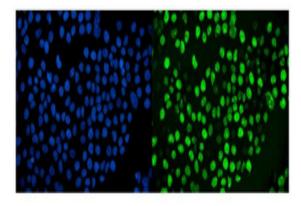


IF analysis of Ki67 using anti-Ki67 antibody. Ki67 was detected in immunocytochemical section of A431 cells. Enzyme antigen retrieval was performed using IHC enzyme antigen retrieval reagent for 15 mins. The cells were blocked with 10% goat serum. And then incubated with 5 μ g/mL mouse anti-Ki67 Antibody overnight at 4°C. DyLight®488 Conjugated Goat Anti-Mouse IgG was used as secondary antibody at 1:100 dilution and incubated for 30 minutes at 37°C. The section was counterstained with DAPI. Visualize using a fluorescence microscope and filter sets appropriate for the label used.


Phone: +1 (415) 906-5211 | Fax: +1 (415) 651-8558






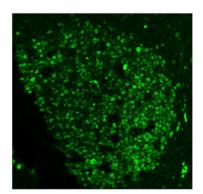

IF analysis of Ki67 using anti-Ki67 antibody. Ki67 was detected in paraffin-embedded section of human tonsil tissue. Heat mediated antigen retrieval was performed in EDTA buffer (pH8.0, epitope retrieval solution). The tissue section was blocked with 10% goat serum. The tissue section was then incubated with 5 µg/mL mouse anti-Ki67 Antibody overnight at 4°C. Biotin conjugated goat anti-mouse IgG was used as secondary antibody and incubated for 30 minutes at 37°C. The tissue section was developed using DyLight®488 Conjugated Avidin. The section was counterstained with DAPI. Visualize using a fluorescence microscope and filter sets appropriate for the label used.

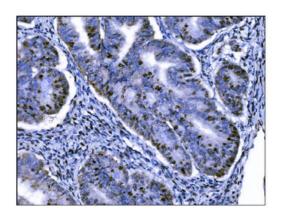

IHC analysis of Ki67 using anti-Ki67 antibody. Ki67 was detected in paraffin-embedded section of human cervical cancer tissue. Heat mediated antigen retrieval was performed in EDTA buffer (pH8.0, epitope retrieval solution). The tissue section was blocked with 10% goat serum. The tissue section was then incubated with 2 μ g/ml mouse anti-Ki67 Antibody overnight at 4°C. Biotinylated goat anti-mouse IgG was used as secondary antibody and incubated for 30 minutes at 37°C. The tissue section was developed using Strepavidin-Biotin-Complex (SABC) with DAB as the chromogen.

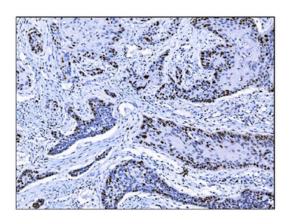


IHC analysis of Ki67 using anti-Ki67 antibody. Ki67 was detected in paraffin-embedded section of human esophageal squamous carcinoma tissue. Heat mediated antigen retrieval was performed in EDTA buffer (pH8.0, epitope retrieval solution). The tissue section was blocked with 10% goat serum. The tissue section was then incubated with 2 μ g/ml mouse anti-Ki67 Antibody overnight at 4°C. Biotinylated goat anti-mouse IgG was used as secondary antibody and incubated for 30 minutes at 37°C. The tissue section was developed using Strepavidin-Biotin-Complex (SABC) with DAB as the chromogen.


Phone: <u>+1 (415) 906-5211</u> | Fax: <u>+1 (415) 651-8558</u>






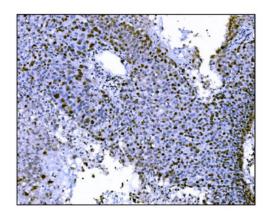

IHC analysis of Ki67 using anti-Ki67 antibody. Ki67 was detected in paraffin-embedded section of human liver tissue. Heat mediated antigen retrieval was performed in EDTA buffer (pH8.0, epitope retrieval solution). The tissue section was blocked with 10% goat serum. The tissue section was then incubated with 2 μ g/ml mouse anti-Ki67 Antibody overnight at 4°C. Biotinylated goat anti-mouse IgG was used as secondary antibody and incubated for 30 minutes at 37°C. The tissue section was developed using Strepavidin-Biotin-Complex (SABC) with DAB as the chromogen.

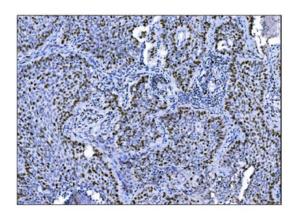

IHC analysis of Ki67 using anti-Ki67 antibody. Ki67 was detected in paraffin-embedded section of human lung cancer tissue. Heat mediated antigen retrieval was performed in EDTA buffer (pH8.0, epitope retrieval solution). The tissue section was blocked with 10% goat serum. The tissue section was then incubated with 2 μ g/ml mouse anti-Ki67 Antibody overnight at 4°C. Biotinylated goat anti-mouse IgG was used as secondary antibody and incubated for 30 minutes at 37°C. The tissue section was developed using Strepavidin-Biotin-Complex (SABC) with DAB as the chromogen.

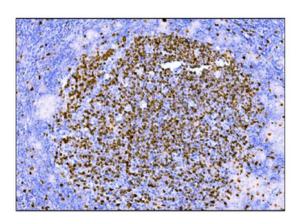

IHC analysis of Ki67 using anti-Ki67 antibody. Ki67 was detected in paraffin-embedded section of human tonsil tissue. Heat mediated antigen retrieval was performed in EDTA buffer (pH8.0, epitope retrieval solution). The tissue section was blocked with 10% goat serum. The tissue section was then incubated with 2 µg/ml mouse anti-Ki67 Antibody overnight at 4°C. Biotinylated goat anti-mouse IgG was used as secondary antibody and incubated for 30 minutes at 37°C. The tissue section was developed using Strepavidin-Biotin-Complex (SABC) with DAB as the chromogen.

Email: $\underline{info@biorbyt.com}$, $\underline{support@biorbyt.com}$ Phone: $\underline{+1 (415) 906-5211}$ | Fax: $\underline{+1 (415) 651-8558}$